

Rabbit anti Bak Polyclonal Antibody Alternative Name(s): BCL2-antagonist/killer 1; BAK, BAK-LIKE, BCL2L7, CDN1

Order Information

- Description: Bak
- Catalogue: 500-3324
- Lot: See label
- Size: 100ug/200ul
- Host: Rabbit
- Clone: nan
- Application: IHC(P)
- Reactivity: Hu, Ms, Hm

ANTIGEN PREPARATION

A synthetic peptide corrsponding to the internal segment of human Bak

BACKGROUND

BAK1 belongs to the BCL2 protein family. BCL2 family members form oligomers or heterodimers and act as anti- or pro-apoptotic regulators that are involved in a wide variety of cellular activities. It localizes to mitochondria, and functions to induce apoptosis. BAK1 interacts with and accelerates the opening of the mitochondrial voltage-dependent anion channel, which leads to a loss in membrane potential and the release of cytochrome c. BAK1 also interacts with the tumor suppressor P53 after exposure to cell stress.

PURIFICATION

The Rabbit IgG is purified by Epitope Affinity Purification

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antbody recognizes human and mice Bak . The other species not tested.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at -20oC to -70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*


- Western Blot: 0.1-1 µg/ml
- ELISA: 0.01-0.1 µg/ml
- Immunoprecipitation: 2-5 µg/ml
- IHC: 2-10 µg/ml
- · Flow cytometry: Not tested
- Molecular Weight: 23.0
- Positive Control: Kidney Tissue
- Cellular Location: Cell Membrane

*Optimal dilutions should be determined by researchers for the specific applications.

FOR RESEARCH USE ONLY.

AbboMax, Inc 2528 Qume Drive, Suite 8, San Jose, California 95131, USA 1 408-573-1898 (Tel). 1 408-573-1858 (Fax). www.abbomax.com info@abbomax.com



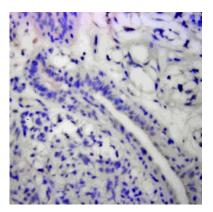

Immunohistochemistry: Human Pancreatic carcinoma (FFPE) stained with Rabbit anti-Bak (Cat# 500-3324) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES