

Rabbit anti V1bR (pY338) Polyclonal Antibody

Alternative Name(s): arginine vasopressin receptor 1B; AVPR1B

Order Information

Description: V1bR (pY338)
Catalogue: 620-530
Lot: See label
Size: 100ug/200ul
Host: Rabbit
Clone: nan

• Application: IHC(P), WB

• Reactivity: Hu

ANTIGEN PREPARATION

A synthetic peptide PWI (pY) MGFN containing a phosphorylation site Tyrosine 338 of human V1bR protein.

BACKGROUND

V1bR, Vasopressin V1b receptor, is a receptor for arginine vasopressin. This receptor belongs to the subfamily of G-protein coupled receptors which includes AVPR1A, V2R and OXT receptors. Its activity is mediated by G proteins which stimulate a phosphatidylinositol-calcium second messenger system. V1bR is primarily located in the anterior pituitary, where it stimulates ACTH release. It is expressed at high levels in ACTH-secreting pituitary adenomas as well as in bronchial carcinoids responsible for the ectopic ACTH syndrome.

PURIFICATION

The Rabbit IgG is purified by site-modified Epitope Affinity Purification.

FORMULATION

This affinity purified antibody is supplied in sterile Tris-buffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes human V1bR (pY338) protein with a phosphorylation site Tyrosine 338.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at -20oC to -70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

• Western Blot: 0.1-1 μg/ml

• ELISA: 0.01-0.1 μg/ml

• Immunoprecipitation: 2-5 µg/ml

• IHC: 2-10 µg/ml

· Flow cytometry: Not tested


• Molecular Weight: 50.0

• Positive Control: Kidney Tissue

Cellular Location: Cell Membrane

^{*}Optimal dilutions should be determined by researchers for the specific applications.



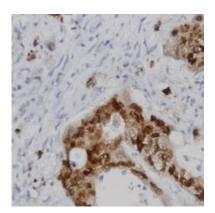

Immunohistochemistry: Human colon carcinoma (FFPE) stained with Rabbit anti-VIR (pY338) (Cat# 520-530) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES