



Rabbit anti BCAR3 (Non-phospho) Polyclonal Antibody

Alternative Name(s): breast cancer anti-estrogen resistance 3

Order Information

- **Description:** BCAR3 (Non-phospho)
- **Catalogue:** 600-090
- **Lot:** See label
- **Size:** 100ug/200ul
- **Host:** Rabbit
- **Clone:** nan
- **Application:** IHC(P), WB
- **Reactivity:** Hu, Rt, Ms

ANTIGEN PREPARATION

A synthetic peptide surrounding of –DRTPE- of human BCAR3 protein WITHOUT any phosphorylation site. This sequence is identical to human, rat, mouse, canis and bovine.

BACKGROUND

A synthetic peptide surrounding of –DRTPE- of human BCAR3 protein WITHOUT any phosphorylation site. This sequence is identical to human, rat, mouse, canis and bovine.

PURIFICATION

The Rabbit IgG is purified by Epitope Affinity Purification

FORMULATION

This affinity purified antibody is supplied in sterile Phosphatebuffered saline (pH7.2) containing antibody stabilizer

SPECIFICITY

This antibody recognizes Non-phosphorylated Threonin 130 of human BCAR3 protein. It cross react to human, mouse and rat.

STORAGE

The antibodies are stable for 24 months from date of receipt when stored at –20oC to –70oC. The antibodies can be stored at 2oC-8oC for three month without detectable loss of activity. Avoid repeated freezing-thawing cycles.

APPLICATIONS/SUGGESTED WORKING DILUTIONS*

- Western Blot: 0.1-1 µg/ml
- ELISA: 0.01-0.1 µg/ml
- Immunoprecipitation: 2-5 µg/ml
- IHC: 2-10 µg/ml
- Flow cytometry: Not tested
- Molecular Weight: 93.0
- Positive Control: Kidney Tissue
- Cellular Location: Cell Membrane


*Optimal dilutions should be determined by researchers for the specific applications.

FOR RESEARCH USE ONLY.

AbboMax, Inc 2528 Qume Drive, Suite 8, San Jose, California 95131, USA
1 408-573-1898 (Tel). 1 408-573-1858 (Fax). www.abbomax.com info@abbomax.com



DATA ATTACHMENTS

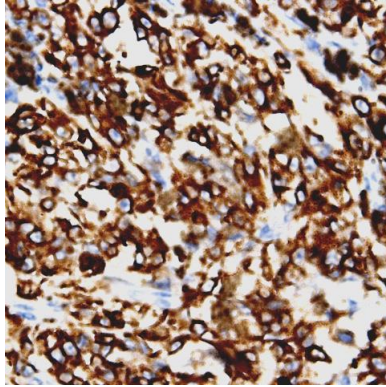

Immunohistochemistry: Human lymph node (FFPE) stained with Rabbit anti-BCAR3 (paired T130) (Cat# 600-090) at 1:200 for 10 min @ RT. Staining of formalin-fixed tissue requires boiling tissue sections in 10 mM Citrate Buffer, pH 6.0 for 10 min followed by cooling at RT for 20 min.

REFERENCES

FOR RESEARCH USE ONLY.

AbboMax, Inc 2528 Qume Drive, Suite 8, San Jose, California 95131, USA
1 408-573-1898 (Tel). 1 408-573-1858 (Fax). www.abbomax.com info@abbomax.com